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#### 2. Characteristics / properties Polyan Plus Economic Line Blanks

- resistant to abrasion
- high tensile strength
- high flexural strength
- plaque resistant and biocompatible
- Color fidelity and high translucency
- Residual monomer content is below 0.5% and free from benzoyl peroxide; therefore suitable for allergy sufferers

#### 3. Composition

Acrylic polymer based on methyl methacrylate

#### 4. Material-mechanical specifications

E-Modul	2800 MPa	DIN EN ISO 527 compliant
Tensile strength	76,3 MPa	DIN EN ISO 527 compliant
Fracture stress	70,3 MPa	DIN EN ISO 527 compliant
Flexural strength	109 MPa	DIN EN ISO 527 compliant
Water absorption	19,36 µg/mm	DIN EN ISO 10477 compliant
Water solubility	< 1µg/mm	DIN EN ISO 10477 compliant

#### 5. Delivery forms of the Blanks

See current CAD/CAM catalog


#### 6. Indication

**Polyan Plus Economic Line** are milling blanks for the production of total prostheses for permanent use.

#### 7. Notes on Design

Minimum wall thickness:

- 0,6 mm (zervikal)
- 1,2 mm (okklusal)
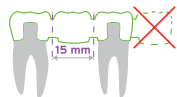
#### 8. Instructions for milling and processing

Milling needs to be done with a milling strategy (forward speed, number of rotation) adapted for the PMMA. Contact the manufacturer of your milling machine for further details. Do not use any diamond-coated tools, which can lead to overheating (luting) of the material.

Milling cutters, which produce good results including adequate water and air colling, are available at Dental Plus. They are called: „Regenbogen-Fräser“ / “Rainbow-Milling-Cutters”. To uncise the milled framework out of the milling-blank a smooth geared carbide milling cutter is used. Separate the connectors carefully and without pressure.

Restorations can be pre-polished with appropriate silicon burnishes and goat-hair-brushes. The mirror finish is carried out with Polirapid-Polishing-Paste.

Scrub / Ø Milling cutter	2 (2,5) mm	Scrub / Ø Milling cutter	1 mm
Rotation (rpm)	19.000 - 22.000	Rotation (rpm)	25.000 - 27.000
Feed rate x/y (mm/min)	1.800 – 2.100	Feed rate x/y (mm/min)	1.200
Feed rate z (mm/min)	1.000	Feed rate z (mm/min)	1.000 - 1.200
Z-lining	0,8 - 1,5 mm	Z-lining	0,05 - 0,08 mm
Scrub / Ø Milling cutter	2 (2,5) mm	Details	0,6 mm
Rotation (rpm)	25.000 - 27.000	Rotation (rpm)	25.000
Feed rate x/y (mm/min)	1.800 - 2.000	Feed rate x/y (mm/min)	2.000 - 2.200
Feed rate z (mm/min)	1.000	Feed rate z (mm/min)	700 - 900
Z-lining	0,2 mm	Z-lining	0,01 - 0,05
Distance to the blank outer edge	2mm	Mandatory milling geometry	1-fluted

#### 9. Customizing

Connecting with resign and composites is possible.

#### 10. Cleaning

High concentrated alcohol and benzene generally lead to cracking (Crackles) in PMMA. Cleaning and disinfection fluids containing these ingredients should not be used.

Steam results in heat and pressure load, which should be avoided. Alternative cleaning possibilities are mechanical methods, or the ultrasound machine.

#### 11. Attachment

For the attachment of restorations made out of **Polyan Plus Economic Line** all proven provisory and permanent cements are suitable.

#### 12. Reactions

Reactions with other products used in dentistry are not to be expected. If reactions in the mouth cannot be excluded, the denture must be removed.

#### 13. Storage

**Polyan Plus Economic Line** Blanks should be stored at room temperature and protected against direct sunlight.

#### 14. Expiration date

Check tag on the packaging.

#### 15. First-Aid Measures

Due to the processing of **Polyan Plus Economic Line** composite Blanks and its properties no health hazards are to be expected. If due to special circumstances any reactions take place, please note the following instructions:

- After inhalation: In case of problems seek medical advice.
- After skin contact: Rinse with water and seek medical advice.
- After eye contact: If mechanical irritation occurs flush eyes with water and seek medical advice.
- After ingestion: In case of discomfort seek medical advice.

#### 16. Warning notes

- **The use of denture finisher can result in embrittlement (breakage) of the material if ingredients are used.**
- **Polyan Plus Economic Line** Blanks can be used in all milling machines. Note, material overheating may be the result of improper use of cutting tools, polishing brushes, steam radiators and water baths. This can lead to damage of the material.